



UK Biomaterials Product catalogue

Valid from January 2024















Content

Introduction	
Biomaterials for hard and soft tissue regeneration	2
Bone graft substitutes	
MinerOss® Blend	3
MinerOss® Cancellous	4
MinerOss® Block Allograft	4
MinerOss® Putty	5
MinerOss® XP	6
MinerOss® X	7
SynMax™	8
Membranes	
Striate+™	10
Mem-Lok® Pliable	11
Mem-Lok® RCM	12
Mem-Lok® Pericardium	13
PermaPro™	14
Cytoplast™ Titanium-Reinforced Dense PTFE Membranes	15 & 16
Wound dressings	
BioPlug and BioStrip	16
Soft Tissue Regeneration	
NovoMatrix®	17
Bone fixation and membrane stabilization	
truFIX System	18
truTACK, truSCREW and truTENT	19
Design de la companya	00.00
Product overview	22–29
Pafarancas	30

Biomaterials for hard and soft tissue regeneration

Bone graft substitutes, membranes, Reconstructive Tissue Matrix and wound dressings

The correct choice of biomaterials is crucial to achieve optimal clinical outcomes - in functional, structural and esthetic terms. Our portfolio of biomaterials offers you a comprehensive range of products for virtually all requirements needed for the regeneration of hard and soft tissue deficits. The product catalogue provides a summary of our entire biomaterials portfolio. It serves as a guide and aid for the selection of suitable biomaterials.

Our product portfolio includes allogeneic (human origin), xenogeneic (porcine and bovine origin) and synthetic bone substitute materials and membranes. Due to their structural properties and manufacturing processes, the materials differ in their resorption behavior as well as their handling.


The range of allogeneic bone substitute material is an allograft made from human donor bone and is subject to high safety standards in the manufacturing process. The range of xenogeneic bone substitute materials is methodically processed from bovine or porcine bone and extensively tested to eliminate potential antigenicity and to provide a favorable environment for new bone growth. Synthetic bone substitute material offers an alternative to commercially available bone substitute materials and extends the treatment spectrum.

In addition to bone substitute materials, our portfolio also includes membranes as well as an acellular dermal tissue matrix of porcine origin. The collagen wound dressings round off our product portfolio allowing complete patient centric planning and treatment selection.







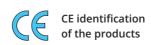





Allograft bone grafting substitute

MinerOss® Blend





In 2004, BioHorizons pioneered the combination of mineralized cortical and cancellous chips. MinerOss Blend is intended to generate an osteoconductive scaffold for volume enhancement while facilitating penetration of blood vessels for revascularization of the site. The Blend now offers a regenerative balance of approximately 70% cortical and 30% cancellous allograft chips.

- mineralized cancellous chips with maximum surface area and early osteoconductive properties1
- mineralized cortical chips for space maintenance and long-term stability¹

applications include:

- extraction site grafting
- periodontal defects
- lateral and vertical ridge augmentation
- sinus grafting
- peri-implant defects
- composite grafting

MinerOss® Blend (Particle size range: 300-1000µm)

Item code	
MIN-C0.5INT	0.5cc Vial
MIN-C1.0INT	1.0cc Vial
MIN-C2.5INT	2.5cc Vial



SEM picture of MinerOss® Blend at 50-fold magnification

Allograft bone grafting substitute

MinerOss® Cancellous





MinerOss Cancellous is a mineralized allograft cancellous bone particle. The osteoconductive properties along with the fast remodeling time allow for rapid revascularization and predictable results.²

MinerOss Cancellous provides ideal handling characteristics for regeneration procedures of all sizes

applications include3:

- ridge and sinus augmentation
- socket grafting

- periodontal defects
- composite grafting

MinerOss® Cancellous (Particle size range: 300-1000µm)

	±.
Item code	
MIN-CAN0.5	0.5cc Vial
MIN-CAN1.0	1.0cc Vial
MIN-CAN2.5	2.5cc Vial

Allograft bone grafting substitute

MinerOss® Block





MinerOss Block Allograft restores bone volume⁴ and provides an alternative to harvesting an autogenous block graft from the patient therefore eliminating the need for a second surgical procedure.

applications include:3

- bone remodeling
- block grafting procedures
- bony defects

MinerOss® Block Alograft

Item code	
MO-BLH10	15 x 10 x 10mm (Illium block)
MO-BLH15	15 x 10 x 15mm (Illium block)

Allograft bone grafting substitute

MinerOss® Putty







MinerOss Putty is uniquely a 100% allograft that includes a powerful combination of mineralized cortical/cancellous chips with demineralized cortical fibres for optimal implant site development.

MinerOss Putty is designed for enhanced regenerative capacity



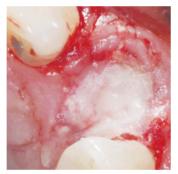



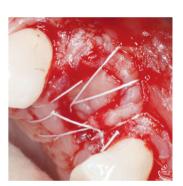

MinerOss Putty case images courtesy of Steve Wallace, DDS, MHS

- over 80% bone graft weight⁵
- 50:50 ratio of mineralized cortical/cancellous chips and demineralized cortical fibers⁵
- all-natural carrier with no extraneous fillers or binding agents
- moldable and cohesive handling characteristics
- time-saving delivery method with reliable clinical outcomes

applications include:3

- ridge and sinus augmentation
- peri-implant defects
- extraction socket grafting
- periodontal defects

MinerOss® Putty (Uniquely 100% allograft)

Item code	
MIN-PT0.5	0.5cc Syringe
MIN-PT1.0	1.0cc Syringe
MIN-PT2.5	2.5cc Syringe
MIN-PT5.0	5.0cc Syringe

Porcine bone graft substitute

MinerOss® XP





MinerOss® XP is a porous bone mineral matrix consisting largely of calcium phosphate. It is obtained by removing organic components from cancellous bone of porcine origin. The inorganic MinerOss® XP bone matrix has macro- and microscopic structures which resemble those of human bone. Due to this trabecular architecture with interconnecting macro- and micropores, the ingrowth of new vessels and bone at the graft site is optimized.

Ideal for following indications

- Augmentation or reconstruction of the alveolar ridge
- Filling of intrabony periodontal defects
- Filling of defects after root resection, apicoectomy, or cystectomy
- Filling of extraction sockets for the protection and preservation of the alveolar ridge
- Sinus floor elevation
- Filling of periodontal defects in conjunction with products for guided tissue regeneration (GTR) or guided bone regeneration (GBR)
- Filling of peri-implant defects in conjunction with products for guided bone regeneration (GBR)

Product features

- Intra and interparticle space ⁶
 - The highly porous structure of MinerOss® XP provides substantial space for the growth of new blood vessels and new bone.
 - More intra and interparticular space is provided for osteoconduction and new bone formation than with comparable
- Rough surface facilitates cell adhesion and spread for bone in-growth6
- High volume fill per unit weight⁶
- Carbonate apatitie substitution promotes better osteoclastic remodelling than hydroxyapatite 7-10

See Product Overview at the end for all the product options and codes for MinerOss® XP



SEM picture of MinerOss® XP at 25-fold magnification macropores andmicropores resemble human bone.

Bovine bone graft substitute

MinerOss® X





MinerOss® X is an anorganic, bovine bone, mineral matrix available in a variety of options. Physically and chemically, the product is comparable to the mineral structure of human bone. The formation and ingrowth of new bone at the implantation site of MinerOss® X is favored, due to its trabecular architecture, interconnecting macro and micro pores and its natural consistency. MinerOss® X Collagen is a combination of 95 % anorganic, cancellous, bovine bone and approximately 5% bovine collagen. This block form allows for convenience during placement and is an ideal solution for many applications, including ridge preservation, minor bone augmentations and periodontal regeneration.

Ideal for following indications

- Augmentation or reconstruction of the alveolar ridge
- Filling of intrabony periodontal defects
- Filling of defects after root resection, apicoectomy, or cystectomy
- Filling of extraction sockets for the protection and preservation of the alveolar ridge
- Sinus floor elevation
- Filling of periodontal defects in conjunction with products for guided tissue regeneration (GTR) or guided bone regeneration (GBR)
- Filling of peri-implant defects in conjunction with products for guided bone regeneration (GBR)

Product features

- Flexible, to meet clinical needs
 - In combination with Mem-Lok® RCM, MinerOss® X preserves ideal space and long-term cell occlusion for maximum bone volume
- Matrix for osseointegration
 - Diffraction patterns are close to the mature native bone diffraction pattern 11
 - High porosity which supports and enhances integration of new bone
- Dependable stability and strength
 - Deproteinized and delipidized, gamma-sterilized
 - Optimal calcium/phosphate balance comparable to human bone 12

See Product Overview at the end for all the product options and codes for MinerOss® X



SEM picture of MinerOss® X at 50-fold magnification macropores and micropores resemble human bone.

Synthetic bone graft substitute

SynMax®





SynMax® is a fully synthetic, safe and biocompatible material that, when brought into an osseous environment, serves as an osteoconductive scaffold to support the ingrowth and fusion of adjacent, vital bone. It's composed of 60 % hydroxyapatite and 40% betatricalcium phosphate.

After implantation the material undergoes a natural remodeling and is gradually resorbed and replaced by new bone. SynMax® is a bone graft material that provides clinicians and their patients with an ideal alternative to human allograft and animal origin bone graft material.^{13–15}

Ideal for following indications

- Sinus lift
- Ridge augmentation
- Intraosseous defects
- Extraction sockets
- Osseous defects
- Furcation defects

SynMax™ (Particle size 500 - 1000µm)


Item code	
BM1013.1005	0.5cc Vial
BM1013.1010	1.0cc Vial

Product features

- 100 % synthetic, no risk of disease transmission, high safety
- Controlled resorption due to biphasic composition
- Very rough surface and high porosity supports integration and bone formation

SynMax[™] (Particle size 800 - 1500 μ m)

Item code	
BM1014.1005	0.5cc Via
BM1014.1020	2.0cc Via

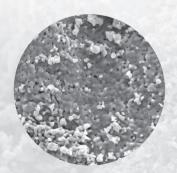

SEM picture of SynMax® at 1000-fold magnification showing microporous structure.



Introducing Xtract Packs: Stronger Together



MinerOss® XP & Striate+™ Xtract Pack	
MinerOss® XP 0.5cc Small & Striate+™ 15 x 20	£130 (£21 saving)
MinerOss® Blend & Striate+™ Xtract Pack	
MinerOss® Blend 0.5cc & Striate+™ 15 x 20	£150 (£25 saving)
MinerOss® Putty & Striate+™ Xtract Pack	
MinerOss® Putty 0.5cc & Striate+™ 15 x 20	£190 (£29 saving)

Introducing Xtract Pack – your solution for individual case management.

BioHorizons Camlog offers the regenerative solutions needed to create a solid base for hard and soft tissue growth – the perfect synergy between natural looking esthetics and long-term function. Now available with a Striate+ $^{\rm IM}$, a noncrosslinked, accelular Type 1 collagen, bi-layered membrane. Our **Xtract Packs** have been designed to simplify the ordering process, case management and stock storage.

Secure your discount now and order an **Xtract Pack** to receive 1 dental membrane and 1 dental bone graft for your next case.

BioHorizons Camlog UK & Ireland

Reflex, Cain road, Bracknell, Berkshire, RG12 8TQ Phone: +44 (0) 1344 752560 • Email: infouk@biohorizons.com www.biohorizonscamlog.com

More information or to enquire:





Introducing Striate+ collagen membrane exceptional clinical results

Non-crosslinked, acellular Type I porcine collagen

does not induce abnormal inflammatory response¹



prevents infiltration of gingival cells while allowing passage of bioactive molecules and proteins¹

Exceptional handling properties

adapts easily to the bone surface, does not collapse or stick when wet

Bi-layer membrane structure

readily conforms to bone surfaces

Bioactive chamber

enables early integration of boneforming cells and provides a favourable environment for osteogenesis¹⁶



CURIOUS?
SCAN THE QR CODE &
REACH OUT TODAY!





Porcine collagen membrane

Striate+™ - natural collagen membrane with a bilayer structure





Striate+™ is a resorbable non-crosslinked collagen barrier membrane for guided bone and tissue regeneration. It is the next generation of collagen membranes which, due to its unique manufacturing process, creates a favorable environment for rapid regeneration of high-quality bone and soft

Optimal for following indications 16

- Treatment of jaw bone defects in oral and maxillofacial surgery
- Guided bone regeneration
- Guided soft tissue regeneration
- Periodontal defect regeneration

Striate+™ Collagen Membrane

Item code	
OCG-152	15 x 20mm
OCG-203	20 x 30mm
OCG-304	30 x 40mm

SEM picture of Striate+ showing the collagen fibre

- Preservation of the collagen architecture of the original tissue due to the optimized manufacturing process.
 - Results in very good processing properties and natural degradation profile
- Exceptional handling properties
 - The membrane can be sutured, screw retained or pinned without tearing or deforming
 - The membrane does not stick in hydrated state
- Bilayer structure to stimulate specific biological responses
 - Osteoconductive, rough side with bioactive chambers serves as a guide for cells and blood vessels
 - faces the bone defect
 - Smooth side with dense collagen structure for barrier function
 - barrier function prevents the migration of gingival cells
 - allows the passage of bioactive molecules
- Minimized inflammatory response
- Predictable results
- Optimal wound stabilization and healing
- Ideal soft tissue integration and hence esthetic results

Porcine collagen membrane

Mem-Lok® Pliable





Mem-Lok® Pliable is a strong and conformable collagen membrane made of highly purified, porcine tissue. Mem-Lok® Pliable offers flexibility and strength. It is easy to handle and simple to fixate. This barrier membrane supports soft tissue and stabilizes the grafting area. Meticulously manufactured from highly purified, intact, porcine collagen and minimally cross-linked, it is biocompatible and predictably resorbable. It is smoothly adaptable to defects and contours and can easily be repositioned. Due to its high suture pullout strength, it can be firmly anchored to the surrounding tissue.

Ideal for following indications

- Augmentation around implants placed in extraction sockets
- Augmentation around implants placed in extended extraction sockets
- Local ridge augmentation for later implantation
- Reconstruction of the alveolar ridge for prosthetic treatment
- Filling of bone defects after root resection, cystectomy, or removal of retained teeth
- Guided bone regeneration in dehiscence defects
- Guided bone regeneration procedures in periodontal defects

Mem-Lok® Pliable Collagen Membrane

Item code	
PBLE-ML1520	15 x 20mm
PBLE-ML2030	20 x 30mm
PBLE-ML3040	30 x 40mm



SEM picture of Mem-Lok® Pliable at 50-fold magnification – not side-specific; dense, uniform single layer

- Special handling characteristics¹⁷
 - Not side-specific
 - Can be placed dry or hydrated
 - Does not adhere to gloves or instruments
 - Simple, easy fixation
 - Single layer, intact collagen
 - Cell occlusive
 - High tear strength
- Supports wound healing¹⁷
 - Reduced degree of inflammation and foreign body response confirmed in pre-clinical testing at early timepoints
 - Protects the graft area from undesirable soft-tissue infiltration uring initial healing phase
 - Enables nutrient transfer
 - Predictable resorption in 12 to 16 weeks
 - Greater initial stability during the critical early weeks of healing due to slow resorption time
- Dependable strength
 - Proven biomechanical strength safeguards fixation
 - In pre-clinical testing, suture pullout strength was three times higher than a comparable collagen membrane.¹⁷

Bovine collagen membrane

Mem-Lok® Resorbable Collagen Membrane (RCM)





Mem-Lok® RCM is manufactured from highly purified, type I bovine collagen. Clinicians can be confident that Mem-Lok® RCM will serve as an effective barrier membrane for bone regeneration. Mem-Lok® RCM supports graft stabilization and bone growth by providing soft tissue support and space maintenance over a predictable timeframe. It is manufactured to ensure predictable resorption rates. Due to its in-vivo stability, it enables easy handling in demanding indications.

Ideal for following indications

- Periodontal defects
- Extraction sockets
- Horizontal ridge enhancement
- Vertical ridge enhancement
- Sinus augmentation
- Dehiscence defects
- Immediate implantation

Mem-Lok® RCM Collagen Membrane

Item code	
RCM-ML1520	15 x 20mm
RCM-ML2030	20 x 30mm
RCM-ML3040	30 x 40mm

SEM picture of Mem-Lok® RCM

- Special handling characteristics 12
 - Membrane only 0.3 mm thick, yet rigid
 - Easy to use due to dimensional stability
 - Easy placement since membrane is not side-specific
 - Potentially reduced treatment time thanks to easy fixation
 - Minimal hydration for optimal bio-adaptability
- Flexible, to meet clinical needs
 - Mem-Lok® RCM maintains ideal space and long-term cell occlusion for maximum bone volume
 - Permeability permits the exchange of essential nutrients during
 - Easily adapts to a whole range of bone defects
- Cell-occlusive for supporting bone regeneration
- Protecting the graft area from undesirable soft tissue infiltration during the initial healing phase
- Predictable resorption after 26 to 38 weeks ¹⁸ eliminates the need of a removal surgery

Dental Membranes

Mem-Lok® Pericardium






Mem-Lok® Pericardium is a long lasting, conformable barrier membrane that drapes easily for graft site contours. It has excellent strength and stability for optimal graft site protection.

applications include: 3


- extraction sockets
- periodontal defects
- block graft coverage
- ridge reconstruction

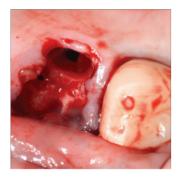



Product features

- tear-resistant, can be sutured 19
- 5-year shelf life⁵
- can be stored at room temperature³
- Predictable resorption time of 16–24 weeks

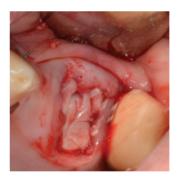



Mem-Lok® Pericardium case images courtesy of Guy Rosenstiel, DMD

Mem-Lok® Pericardium

Item code	
PERI-ML1520INT	15 x 20mm
PERI-ML2030INT	20 x 30mm
PERI-ML3040INT	30 x 40mm

Synthetic PTFE membrane

PermaPro®





PermaPro® is an exceptionally thin, non-resorbable, temporary implantable and biocompatible membrane. It is composed of biologically inert, highdensity polytetrafluoroethylene (PTFE), which acts as an efficient barrier against bacterial and cellular penetration, and can therefore be used for open healing in certain indications.


Ideal for following indications

- For the regeneration of extraction sockets (Socket and Ridge Preservation)
- For use as a space-creating barrier in guided bone regeneration (GBR) and guided tissue regeneration (GTR)
- For covering bone defects during surgical procedures in periodontology, oral and maxillofacial surgery, oral surgery, and implant dentistry

PermaPro® (synthetic PTFE membrane)

Item code	
BM.2005.1520	15 x 20mm
BM.2005.2030	20 x 30mm
BM.2005.3040	30 x 40mm

- 100 % synthetic PTFE barrier membrane
- Ultra-thin (approx. 0.08mm)
- Impervious to bacteria due to dense structure
- Easily removable due to minimal tissue ingrowth into the surface structure
- No need for primary soft tissue closure (indication-dependent) ^{20, 21}
- Easy recovery thanks to blue color
- Rounded edges for minimal tissue trauma
- Easy fixation with sutures or pins
- Higher dimensional stability compared to commercially available collagen membranes
- Augmentation outside the ridge contour
- Synthetic nature no religious or dietary conflicts
- Exposure situations where primary wound closure is not desired (indication dependent)

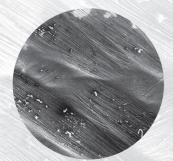

SEM picture of PermaPro® at 30-fold magnification

Dental Membranes

Cytoplast™ Titanium-Reinforced Dense PTFE Membranes





The traditional frame design, incorporating delicate and strategicallyplaced titanium "struts", has more than 25 years of clinical history and successful use in guided bone regeneration. Cytoplast Ti-250 membranes provide a wide range of coverage solutions for cases involving extraction sites, bony defects and ridge augmentation (250 microns thick).

- non-resorbable material allows the clinician to dictate healing time
- passive fit with no memory retention
- lightweight framework is easy to trim and is compliant with overlying soft tissues

Cytoplast™ Titanium-Reinforced Dense PTFE Membranes





OG-TI250AS-2 Ti-250 Anterior Singles (pack of 2) 14mm x 24mm Coverage of single-tooth extraction sites, especially where one or more bony walls are missing.



OG-TI250ATC-2 Ti-250 Anterior Trans Crestal (pack of 2) 24mm x 38mm Designed for bony defects between adjacent teeth, including ridge augmentation.



OG-TI250BL-2	Ti-250 Buccal (pack of 2)	17mm x 25mm
Treatment of large buccal defects.		



OG-TI250PD-2	Ti-250 Posterior Distal (pack of 2)	38mm x 38mm
Designed for large bony defects, including distal extension of the posterior ridge.		



OG-TI250PL-2	Ti-250 Posterior Large (pack of 2)	25mm x 30mm
Suited for treating large b	oony defects, including ridge augmentation.	



OG-TI250PS-2	Ti-250 Posterior Singles (pack of 2)	20mm x 25mm
Suited for covering posterior extraction sites and limited ridge augmentation.		



OG-TI250PST-2	Ti-250 Posterior Singles T2 (pack of 2)	25mm x 36mm
Designed for grafting extraction sites and limited ridge augmentation.		

Dental Membranes

Cytoplast™ Titanium-Reinforced Dense PTFE Membranes





Cytoplast™ Titanium-Reinforced Dense PTFE Membranes

Item Code	Description	Size
OG-TI250PTC-2 Ti-250 Posterior Trans Crestal (pack of 2)		38mm x 38mm
Designed for large bony defects between adjacent teeth, including ridge augmentation.		



OG-TI250XL-2	Ti-250 XL (pack of 2)	30mm x 40mm
Sized to cover very large bony defects, including ridge augmentation. Larger titanium spans more of the PTFE membrane for additional rigidity.		arger titanium spans



OG-TI250XLK-2	Ti-250 XLK (pack of 2)	30mm x 40mm
Sized to cover large bony allows for greater versati	defects, including ridge augmentation. Smalle lity when shaping.	r titanium frame

Bovine collagen wound dressings

BioPlug and BioStrip





BioPlug (Resorbable collagen plug)

Item code BIOPLUG 1cm x 2cm

BioStrip (Resorbable collagen strip)

Item code **BIOSTRIP** 2.5cm x 7.5cm

BioPlug and BioStrip are wound dressings made from bovine collagen. They are designed to absorb blood or fluids and to protect the wound, thus supporting optimal healing. Collagen supports the formation of the blood coagulum and contributes to a rapid stabilization of the wound area.²² Because of their haemostyptic effect, collagen wound dressings are used for the stabilization of extraction sockets and biopsy sampling points as well as in the treatment of smaller wounds.

applications include

- extraction sockets
- denture sores
- oral ulcers

- burns
- surgical wounds
- traumatic wounds

- Fully resorbable in 10 to 14 days
- 10 units per pack
- Packaged sterile

Soft Tissue Regeneration

NovoMatrix[™] Reconstructive Tissue Matrix





NovoMatrix[™] is an acellular dermal matrix derived from porcine tissue. In surgical application, the tear-resistant and easy-to-handle ^{23,24} matrix is an excellent alternative to autologous connective tissue grafts (CTG). There is no need for an intraoral surgical donor site, which reduces morbidity for the patient.

Owing to the manufacturing process, the matrix is free of donor cells. At the same time, the structure of the source tissue remains virtually unchanged, thus supporting the ingrowth of cells and micro-vessels. Proprietary tissue processing enables optimal cell repopulation and revascularization through gentle preparation, resulting in esthetic soft tissue regeneration ²⁹. NovoMatrix™ is supplied pre-hydrated in a patented aqueous phosphatebuffered solution containing matrix stabilizers and can therefore be used promptly without requiring extensive rehydration ³.

Ideal for following indications³

- Increase in attached tissue around teeth and implants
- Reconstruction of the alveolar ridge for prosthetic restoration
- Guided tissue regeneration in recession defects for root coverage

Product features

- The LifeCell™ tissue preparation process results in rapid revascularization.
- Consistent tissue thickness at all times
- Pre-hydrated ready-to-use out of the package following a 2-minute soak in sterile saline or lactated Ringer's solution³
- Store at -8 °C to +30 °C ³

NovoMatrix™ Reconstructive Tissue Matrix

Item code	200
NOVO-1.5X1.5	1.5 x 1.5cm
NOVO-1.5X2.5	1.5 x 2.5cm
NOVO-1.5X4.5	1.5 x 4.5cm

Advantages of NovoMatrix[™] application

Shorter surgery time

The ready-to-use collagen matrix shortens surgery time by eliminating the need for a second donor site ²⁶.

Lower patient morbidity

Avoiding a donor site on the palate eliminates the post-operative pain associated with a second procedure ^{26–28}.

Excellent tissue integration

The application of NovoMatrix[™] supports rapid revascularization, cellular repopulation and minimal inflammatory reactions ^{25, 29–31}.

Natural tissue and color structure

The application of NovoMatrix™ demonstrates irritation-free healing and very good adaptation of the color and tissue structure to the natural surrounding tissue 32.

Rapid and complication-free healing of soft tissue

The application of NovoMatrix[™] supports a positive immunological reaction as well as tissue integration and regeneration 25, 30, 31, 33.



Further information, videos and clinical case studies at www.biohorizonscamlog.com/novomatrix



Bone fixation and membrane stabilization

truFIX System



The truFIX System, distributed by BioHorizons Camlog, offers you everything you need for the fixation of bone blocks and plates as well as the stabilization of membranes. This system incorporates all the necessary components to pick up and drive the truSCREW and truTACK. The truFIX system eliminates the need for multiple systems and unnecessary components, making it user-friendly for your practice.

Product features

- True centered locking connection
- True axial alignment with pickup each time
- truSCREW: self-drilling screw developed for easy insertion with maximum fixation
- truSCREW: Removal Sleeve for disengaging the screw from the driver without damaging screw head
- truTACK: hexagon driver head and barbed tip to pierce without drilling for simple tack insertion and easy screw-like removal
- Easy insertion and easy screw-like removal of the truTACK

The truFIX System includes:

- 1 truFIX Tray (empty)
- 2 truFIX Driver Handle, 98 mm (3.875") long
- 3 truTACK Driver Tip (includes blue Tip Cover)
- 4 truSCREW Driver Tip
- 5 truSCREW Driver Removal Sleeve
- 6 truSCREW Driver Tip, Contra-Angle
- 7 CA 2-Step Countersink Bur, 0.8 mm and 1.6 mm Steps
- 8 1.1 mm Pilot Twist Drill, 29.8 mm
- 9 Pilot Bur, 0.45 mm, 27 mm long, Contra Angle
- 10 External Hex Hand Driver, 0.88 mm
- 11 truSCREW Packaging Removal Tool

(also available separately)

Optional:

truFIX Small Driver Handle, 89 mm (3.5") long



See Product Overview at the end for all the product options and codes for TruFIX.

truTACK, truSCREW and truTENT



The truTACK makes the stabilization of membranes quick and troublefree. Our unique tack incorporates a hexagon on its head and threads on its shaft, allowing for easy removal. The truTACK is placed like a tack and removed like a screw, a feature that you will not find in any other system.

The truSCREW, with its aggressive cutting flutes, is the ideal bone screw for the fixation of small bone within the oral and maxillofacial environment. These cutting flutes (in most instances) eliminate the need for any pre-drilling. The design of the screw ensures an effortless insertion into all types of bone.

The truTENT is a refinement of the truSCREW. Its raised collar and wider head is designed to support a membrane or titanium mesh during augmentation procedures.

truTACK - Bone Tacks (head Ø 2.5 mm)

- Thread Ø 0.7 mm / Total length 3.0 mm (pack of 10)
- Thread Ø 0.7 mm / Total length 5.0 mm (pack of 10)

truSCREW - Bone Screws (head Ø 3.0 mm)


- Thread Ø 1.2 mm / Total length 4.5 mm (pack of 5)
- Thread Ø 1.2 mm / Total length 6.0 mm (pack of 5)
- Thread Ø 1.2 mm / Total length 7.5 mm (pack of 5)
- Thread Ø 1.2 mm / Total length 9.0 mm (pack of 5)
- Thread Ø 1.2 mm / Total length 10.5 mm (pack of 5)
- Thread Ø 1.5 mm / Total length 6.0 mm (pack of 5)
- Thread Ø 1.5 mm / Total length 7.5 mm (pack of 5)
- Thread Ø 1.5 mm / Total length 9.0 mm (pack of 5)
- Thread Ø 1.5 mm / Total length 10.5 mm (pack of 5)
- Thread Ø 1.5 mm / Total length 12.0 mm (pack of 5)
- Thread Ø 1.5 mm / Total length 13.5 mm (pack of 5)
- Thread Ø 1.5 mm / Total length 15.0 mm (pack of 5)
- Thread Ø 2.0 mm / Total length 6.0 mm (pack of 5)
- Thread Ø 2.0 mm / Total length 7.5 mm (pack of 5)
- Thread Ø 2.0 mm / Total length 9.0 mm (pack of 5)
- Thread Ø 2.0 mm / Total length 10.5 mm (pack of 5)
- Thread Ø 2.0 mm / Total length 12.0 mm (pack of 5) Thread Ø 2.0 mm / Total length 13.5 mm (pack of 5)
- Thread Ø 2.0 mm / Total length 15.0 mm (pack of 5)
- truTENT Tenting Screws (head Ø 5.0 mm)

- Thread Ø 1.5 / Total length 10.0 mm / Collar height 4.0 mm
- Thread Ø 1.5 / Total length 12.0 mm / Collar height 6.0 mm
- Thread Ø 1.5 / Total length 14.0 mm / Collar height 8.0 mm

Titanium Meshes





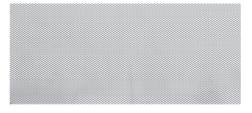







For the reconstruction of extensive combined bony alveolar ridge defects, the use of titanium meshes is advantageous. They serve as a cage to preserve the space created for the augmentate for regeneration. The meshes are adapted intraoperatively to the defect, filled with augmentation material and fixed in positionally stable with screws. They have no barrier function. The titanium meshes are available in different sizes and structures as flat meshes. Depending on the indication, implantation can be performed on one or two sides.

- Titanium Micro Mesh, 120 × 60 mm, 0.1 mm thick
- Titanium Micro Mesh, 34 × 25 mm, 0.1 mm thick
- Titanium Micro Mesh, 152 × 66 mm, 0.2 mm thick
- Titanium Single Butterfly Tenting Mesh, 30 × 80 mm, 0.25 mm thick
- Titanium Tenting Mesh, 13 × 33 mm, 0.2 mm thick
- Titanium Custom Tenting Mesh, 13 × 33 mm, 0.2 mm thick
- Titanium Tenting Mesh, 7 × 14 mm, 0.2 mm thick

See Product Overview at the end for all the product options and codes for TruTACK, TruSCREW, TruTENT and Titanium Meshes.





Bone graft substitutes



Allograft

MinerOss® Blend (Particle size range: 300-1000µm)

Item Code	Volume
MIN-C0.5INT	0.5cc Jar
MIN-C1.0INT	1.0cc Jar
MIN-C2.5INT	2.5cc Jar

MinerOss® Cancellous (Particle size range: 300-1000µm)

Item Code	Volume	
MIN-CAN0.5	0.5cc Vial	
MIN-CAN01.0	1.0cc Vial	
MIN-CAN02.5	2.5cc Vial	

MinerOss® Block Allograft

Item Code	Product size
MO-BLH10	15 x 10 x 10 mm
MO-BLH15	15 x 10 x 15 mm

MinerOss® Putty (100% pure allograft)

Item Code	Volume
MIN-PT0.5	0.5cc Syringe
MIN-PT1.0	1.0cc Syringe
MIN-PT2.5	2.5cc Syringe
MIN-PT5.0	5.0cc Syringe



SynMax® (synthetic bone graft substitute)

Item Code	Volume	Particle size
BM1013.1005	0.5 cm ³	500-1000 μm
BM1013.1010	1.0 cm ³	500–1000 μm
BM1014.1005	0.5 cm ³	800–1500 μm
BM1014.1020	2.0 cm ³	800-1500 μm











Bone graft substitutes



MinerOss® XP Cancellous Syringe (Applicator)

Item Code	Volume	Particle size
MINXP-SYR0.5	0.5 cm ³	250–1000 μm

MinerOss® XP Cancellous (porcine bone graft substitute)

Item Code	Volume	Particle size
MINXP-CAN0.5SM	0.5cc	250-1000 μm
MINXP-CAN1.0SM	1cc	250-1000 μm
MINXP-CAN2.0SM	2cc	250-1000 μm
MINXP-CAN4.0SM	4cc	250-1000 μm
MINXP-CAN1.0LG	1cc	1000-2000 μm
MINXP-CAN2.0LG	2cc	1000–2000 μm





MinerOss® X Cancellous (bovine bone graft substitute)

Item Code	Weight / Volume	Particle size
MINX-CAN0.25GR	0.25 g / 0.6 cm ³	250-1000 µm
MINX-CAN0.5GR	0.5 g / 1.2 cm ³	250-1000 µm
MINX-CAN1.0GR	1.0 g / 2.4 cm ³	250-1000 µm
MINX-CAN2.0GR	2.0 g / 4.7 cm ³	250-1000 µm
MINX-CAN0.25GRL	0.25 g / 0.9 cm ³	1000-2000 μm
MINX-CAN0.5GRL	0.5 g / 1.7 cm ³	1000-2000 μm
MINX-CAN1.0GRL	1.0 g / 3.4 cm ³	1000-2000 μm
MINX-CAN2.0GRL	2.0 g / 6.8 cm ³	1000-2000 μm



MinerOss® X Cancellous Syringe (Applicator)

Item Code	Volume	Particle size
MINX-SYR0.5	0.5 cm ³	250–1000 μm

Membranes



Striate+[™] (porcine collagen membrane)

Item Code	Product size
OCG-152	15 × 20 mm
OCG-203	20 × 30 mm
OCG-304	30 × 40 mm



Mem-Lok® Pliable (porcine collagen membrane)

Item Code	Product size
PBLE-ML1520	15 × 20 mm
PBLE-ML2030	20 × 30 mm
PBLE-ML3040	30 × 40 mm





Mem-Lok® RCM (bovine collagen membrane)

Item Code	Product size
RCM-ML1520	15 × 20 mm
RCM-ML2030	20 × 30 mm
RCM-ML3040	30 × 40 mm





Allograft

Mem-Lok Pericardium (human collagen membrane)

Item Code	Product size
PERI-ML1520INT	15 × 20 mm
PERI-ML2030INT	20 × 30 mm
PERI-ML3040INT	30 × 40 mm



15 × 20 mm

20 × 30 mm

30 × 40 mm

Membranes



Synthetic

PermaPro® (synthetic PTFE membrane)

Art. No.	Product size
BM2005.1520	15 × 20 mm
BM2005.2030	20 × 30 mm
BM2005.3040	30 × 40 mm



15 × 20 mm



Cytoplast™ Titanium-Reinforced Dense PTFE Membranes

Item Code	Item Description	Product size
OG-TI250ANL-2	Ti-250 Anterior Narrow (pack of 2)	12mm x 24mm
Coverage of narrow single-tooth extraction sites, especially where one bony wall is missing.		

OG-TI250AS-2	Ti-250 Anterior Singles (pack of 2)	14mm x 24mm
Coverage of single-tooth extraction sites, especially where one or more bony walls are missing.		more bony walls

OG-TI250ATC-2	Ti-250 Anterior Trans Crestal (pack of 2)	24mm x 38mm
Designed for bony defe	cts between adjacent teeth, including ridge	augmentation.

OG-TI250BL-2	Ti-250 Buccal (pack of 2)	17mm x 25mm
Treatment of large buc	cal defects.	

OG-TI250PD-2	Ti-250 Posterior Distal (pack of 2)	38mm x 38mm
Designed for large bony defects, including distal extension of the posterior ridge.		

OG-TI250PL-2	Ti-250 Posterior Large (pack of 2)	25mm x 30mm
Suited for treating large	bony defects, including ridge augmentation	n.

OG-TI250PS-2	Ti-250 Posterior Singles (pack of 2)	20mm x 25mm
Suited for covering posterior extraction sites and limited ridge augmentation.		

OG-TI250PST-2	Ti-250 Posterior Singles T2 (pack of 2)	25mm x 36mm
Designed for grafting ex	traction sites and limited ridge augmentation	

OG-TI250PTC-2	Ti-250 Posterior Trans Crestal (pack of 2)	38mm x 38mm
Designed for large hopy defects between adjacent teeth, including ridge augmentation		

OG-TI250XL-2	Ti-250 XL (pack of 2)	30mm x 40mm
Sized to cover very large bony defects, including ridge augmentation. Larger titanium		

OG-TI250XLK-2	Ti-250 XLK (pack of 2)	30mm x 40mm	
	Sized to cover large bony defects, including ridge augmentation. Smaller titanium frame allows for greater versatility when shaping.		



Membranes



Synthetic





Ti-250 Anterior Narrow 12mm x 24mm



Ti-250 Anterior Singles 14mm x 24mm



Ti-250 Anterior Trans Crestal 24mm x 38mm



Ti-250 Buccal 17mm x 25mm




Ti-250 Posterior Distal 38mm x 38mm



Ti-250 Posterior Large 25mm x 30mm



Ti-250 Posterior Singles 20mm x 25mm

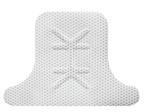

Ti-250 Posterior Singles T2 25mm x 36mm



Ti-250 Posterior Trans Crestal 38mm x 38mm



Ti-250 XL 30mm x 40mm



Ti-250 XLK 30mm x 40mm

Cytoplast Titanium-Reinforced Dense PTFE Membrane size reference

Wound dressings



BioPlug and BioStrip (bovine collagen wound dressings)

Item Code	Product size	
BIOPLUG	10 × 20 mm	Pack
BIOSTRIP	25 × 75 mm	Pack

of 10 of 10





 25×75 mm (BioStrip size reference)



1cm x 2cm (BioPlug size)

Product overview

Soft Tissue Regeneration



NovoMatrix® (porcine, acellular dermal matrix)

Item Code	Product size
NOV1515	15 × 15 mm
NOV1525	15 × 25 mm
NOV1545	15 × 45 mm
NOV2545	25 × 45 mm



15 × 25 mm

15 × 45 mm



25 × 45 mm

Bone fixation and membrane stabilization

truFIX System (starter set)

Item Code	Product Description
45418015K	truFIX System (consisting of Tray / Instruments 1–11) truTACK 3.0mm (10 pack) truSCREW 1.2mm / 9mm (5 pack) truSCREW 1.2mm / 10.5mm (5 pack) truSCREW 1.5mm / 9mm (5 pack)



truFIX System (complete)

Item Code	Product Description
45418015	truFIX System (consisting of Tray / Instruments 1–11)

Tray / Instruments truFIX System

Item Code	Product Description
45418501	1 truFIX Tray (empty)
45419001	2 truFIX Driver Handle, 98 mm (3.875") long
45417001	3 truTACK Driver Tip (includes blue Tip Cover)
45415001	4 truSCREW Driver Tip
45417901	5 truSCREW Driver Removal Sleeve
45415201	6 truSCREW Driver Tip, Contra-Angle
45440203	7 CA 2-Step Countersink Bur, 0.8 mm and 1.6 mm Steps
45440202	8 1.1 mm Pilot Twist Drill, 29.8 mm
502700045	9 Pilot Bur, 0.45 mm, 27 mm long, Contra Angle
20157702	10 External Hex Hand Driver, 0.88 mm
45450201	11 truSCREW Packaging Removal Tool
45419501	truFIX Small Driver Handle, 89 mm (3.5") long (optional)



Titanium Meshes (Grade 1, non-sterile, non-resorbable)

Item Code	Product Description
00039429	Titanium Micro Mesh, 34 × 25 mm, 0.1 mm thick
00039430	Titanium Micro Mesh, 120 × 60 mm, 0.1 mm thick
00039433	Titanium Micro Mesh, 152 × 66 mm, 0.2 mm thick
00039434	Titanium Tenting Mesh, 7 \times 14 mm, 0.2 mm thick
00039440	Titanium Tenting Mesh, 13 × 33 mm, 0.2 mm thick
00039442	Titanium Custom Tenting Mesh, 13 × 33 mm, 0.2 mm thick
00039444	Titanium Single Butterfly Tenting Mesh, 30×80 mm, 0.25 mm thick



Bone fixation and membrane stabilization

truTACK (Bone Tacks, head Ø 2.5 mm, sterile, single-use)

Item Code	Thread Ø / Total length	Pack size
09600313	0.7 mm / 3.0 mm	Pack of 10
09600314	0.7 mm / 5.0 mm	Pack of 10



truSCREW (Bone Screws, head Ø 3.0 mm, sterile, single-use)

Item Code	Thread Ø / Total length	Pack size
45427202	1.2 mm / 4.5 mm	Pack of 5
45427203	1.2 mm / 6.0 mm	Pack of 5
45427204	1.2 mm / 7.5 mm	Pack of 5
45427205	1.2 mm / 9.0 mm	Pack of 5
45427206	1.2 mm / 10.5 mm	Pack of 5
45427502	1.5 mm / 6.0 mm	Pack of 5
45427503	1.5 mm / 7.5 mm	Pack of 5
45427504	1.5 mm / 9.0 mm	Pack of 5
45427505	1.5 mm / 10.5 mm	Pack of 5
45427506	1.5 mm / 12.0 mm	Pack of 5
4542.7507	1.5 mm / 13.5 mm	Pack of 5
45427508	1.5 mm / 15.0 mm	Pack of 5
45428002	2.0 mm / 6.0 mm	Pack of 5
45428003	2.0 mm / 7.5 mm	Pack of 5
45428004	2.0 mm / 9.0 mm	Pack of 5
45428005	2.0 mm / 10.5 mm	Pack of 5
4542.8006	2.0 mm / 12.0 mm	Pack of 5
45428007	2.0 mm / 13.5 mm	Pack of 5
45428008	2.0 mm / 15.0 mm	Pack of 5



truTENT (Tenting Screws, head Ø 5.0 mm, sterile, single-use)

Item Code	Thread Ø / Total length	Collar height
45431001	1.5 mm / 10.0 mm	4.0 mm
45431002	1.5 mm / 12.0 mm	6.0 mm
45431003	1.5 mm / 14.0 mm	8.0 mm



References

- [1] The Efficacy of Mineralized Allograft Cortical and Cancellous Chips in Maxillary Sinus Augmentations. Nevins M, Parma Benfenati S, Janke U, Kleyer A, Rasperini G, Tinti C, Schupbach P, Kim D. Int | Periodontics Restorative Dent 2014;34:789-793
- [2] The efficacy of mineralized allograft cortical and cancellous chips in maxillary sinus augmentations. Nevins M, Janke U, Rasperini G, Schupbach P. Int | Periodontics Restorative Dent.
- [3] Reference manufacturer's Instructions for Use (IFU) package insert.
- Comparison of bone resorption rates after intraoral block bone and guided bone regeneration augmentation for the reconstruction of horizontally deficient maxillary alveolar ridges. Gultekin B. Bedeloglu E. Kose T. Mijiritsky E. BioMed Research International. 2016: Article ID 4987437.
- [6] Data on file, Shu-Thung Li, Ph. D. et al.: Isolation and Characterization of a Porous Carbonate Apatite From Porcine Cancellous Bone. Science, Technology, Innovation, Aug. 2014: 1–13.
- [7] Spense G, Patel N, Brooks R, Rushton N: Osteoclastogenesis on hydroxyapatite ceramics: the effect of carbonate substitution. J Biomed Mater Res A., Mar 15, 2010; 92(4):1292–300.
- [8] Ellies LG, Carter JM, Natiella JR, Featherstone JDB, Nelson DGA: Quantitative Analysis of Early In Vivo Tissue Response to Synthetic Apatite Implants. | Biomed Mater, 1988, Res 22:137–148.
- [9] Landi E, Celotti G, Logroscino G, Tampieri A: Carbonated Hydroxyapatite as Bone Substitute. Journal of the European Ceramic Society, 2003, 23:2931–2937.
- [10] Spense G, Patel N, Brooks R, Rushton N: Carbonate Substituted Hydroxyapatite: Resorption by Osteoclasts Modifies the Osteoblastic Response. Journal of Biomedical Materials Research, 2009, Part A 217-224
- [11] Shu-Tung Li, Hui-Chen Chen and Debbie Yuen: Comparison of a New Natural Bovine Bone Mineral (Carbonate Apatite Anorganic Bone) to Currently Marketed NuOss™ and Bio-Oss®: In Vitro and In Vivo Evaluations. Collagen Matrix, Inc., Oakland, New Jersey 07436.
- [12] Gonshor A. Chris L Tye: Evaluation of Anorganic Boyine Bone Mineral in Post-extraction Alveolar Sockets: A Case Series, Journal of Osseointegration, March 2010; 1(2).
- [13] Binderman et al. Tissue Engineering of Bone: Critical Evaluation of Scaffold Selection. Haim Tal, IntechOpen. April 4th 2012.
- [14] Jelusic et al. Monophasic ß-TCP vs. biphasic HA/ß-TCP in two-stage sinus floor augmentation procedures a prospective randomized clinical trial. Clin Oral Implants Res. 2017 Oc t:28(10):e175-e183.
- [15] Lorenz et al. Investigation of peri-implant tissue conditions and peri-implant tissue stability in implants placed with simultaneous augmentation procedure: a 3-year retrospective followup analysis of a newly developed bone level implant system. Int J Implant Dent. 2017 Sep 5;3(1):41.
- [16] Allan B. et al. Collagen Membrane for Guided Bone Regeneration in Dental and Orthopedic Applications. Tissue Engineering 2020.
- [17] Data on file, Li ST, Yuen D, Martin D, Lee NS: A comparative study of a new porcine collagen membrane to BioGide®. Science, Technology, Innovation. February 1–5, 2015.
- [18] Data on file, Debbie Yuen et al.: Prediction of in vivo stability of a resorbable, reconstituted type I collagen membrane by in vitro methods. World Biomaterials Congress Transactions, Sixth World Biomaterials Congress Transactions. Collagen Matrix Inc., Franklin Lakes, NJ 07417 USA
- [19] Mechanical testing of pericardium for manufactring prosthetic heart valvles. Aguiari P et al. Interactive CardioVascular and Thoracic Surgery. 22 (2016) 72-84.
- [20] Zafiropoulos et al. Open-Healing Socket Preservation with a Novel Dense Polytetrafluoroethylene (dPTFE) Membrane: A Retrospective Clinical Study. Medicina (Kaunas). 2020 Apr
- [21] Papi et al. The Use of a Non-Absorbable Membrane as an Occlusive Barrier for Alveolar Ridge Preservation: A One Year Follow-Up Prospective Cohort Study. Antibiotics (Basel). 2020 Mar 3;9(3):110.
- [22] Nuvttens BP et al. Platelet adhesion to collagen, Thromb Res. 2011 Jan: 127.
- [23] Data on file, Allergan. NovoMatrix™ Mechanical testing, Preclinical Data
- [24] Data on file, Allergan, INT/0204/2018.
- [25] Suárez-López Del Amo F, Rodriguez JC, Asa'ad F, Wang HL. Comparison of two soft tissue substitutes for the treatment of gingival recession defects: an animal histological study. | Appl Oral Sci., 2019;27:e20180584.
- [26] Griffin T, Cheung W, Athanasios Z, Damoulis P. Postoperative Complications Following Gingival Augmentation Procedures. J Periodontology 2006;77:2070-2079.
- [27] Aguirre-Zorzano LA, García-De La Fuente AM, Estefanía-Fresco R, Marichalar-Mendía X. Complications of harvesting a connective tissue graft from the palate. A retrospective study and description of a new technique. J Clin Exp Dent. 2017;9(12):e1439-45.
- [28] Tavelli L, Asa'ad F, Acunzo R, Pagni G, Consonni D, Rasperini G. Minimizing Patient Morbidity Following Palatal Gingival Harvesting: A Randomized Controlled Clinical Study. The International Journal of Periodontics & Restorative Dentistry 38(6):e127-e134 \cdot November 2018.
- [29] Harper JR, McQuillan DJ. Extracellular wound matrices: a novel regenerative tissue matrix (RTM) technology for connective tissue reconstruction. Wounds. 2007;19(6):163-168.
- [30] Sandor M, Leamy P, Assan P, et al. Relevant in vitro predictors of human acellular dermal matrix-associated inflammation and capsule formation in a nonhuman primate subcutaneous tissue expander model. Eplasty. 2017;17:e1-e21.
- [31] Xu H, Wan H, Sandor M, et al. Host response to human acellular dermal matrix transplantation in a primate model abdominal wall repair. Tissue Eng Part A. 2008;14(2):2009-2019.
- [32] Van Orten A. Peri-implant thickening of soft tissue stable and functional. Implantologie Journal 5 | 2020.
- [33] Sandor M, Xu H, Connor J, et al. Host response to implanted porcine-derived biologic materials in a primate model of abdominal wall repair. Tissue Eng Part A. 2008;14(12):2021-2031.

Shop online 24/7

convenient shopping on your schedule

Explore our online store and experience the easy to use customer interface and benefit from:


Easy and efficient, 24/7 ordering
FREE Next Day Shipping*
Order tracking
Inventory availability
Visibility of invoices and statements
Simple reordering
Back order reporting
eReturns portal

Visit

store.biohorizons.com

and register today!



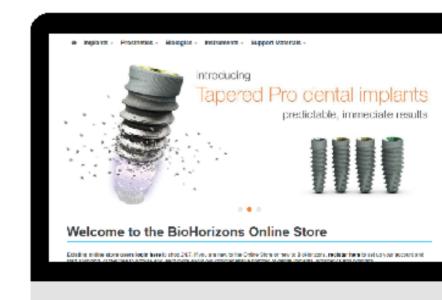



BioHorizons Camlog | Reflex, Cain Road | Bracknell, Berkshire | RG12 1HL | United Kingdom infouk@biohorizons.com | www.biohorizonscamlog.com

Customer Service

Phone +44 (0) 1344 752560 | Email ordersuk@biohorizons.com | Shop store.biohorizons.com

Headquarters

CAMLOG Biotechnologies GmbH | Margarethenstr. 38 | 4053 Basel | Switzerland Phone +41 61 565 41 00 | Fax +41 61 565 41 01 | info@camlog.com | www.biohorizonscamlog.com

MinerOss® X, MinerOss® XP, Mem-Lok® RCM, Mem-Lok® Pliable, BioPlug and BioStrip are manufactured by Collagen Matrix, Inc. MinerOss® Putty, MinerOss® Cancellous, MinerOss® Block, and Mem-Lok® Pericardium are processed by DCI Donor Services Tissue Bank. MinerOss® Cortical & Cancellous Chips are processed by ACE Surgical Supply. SynMax® and PermaPro® are manufactured by botiss biomaterials GmbH. NovoMatrix® is manufactured by LifeCell™ Corporation, an Allergan affiliate. truFIX, truTACK, truSCREW, truTENT and the Titanium Meshes are manufactured by ACE Surgical Supply Co., Inc.

BioHorizons®, MinerOss®, Mem-Lok® and NovoMatrix® are registered trademarks of BioHorizons. Cytoplast® is a registered trademark of Osteogenics Biomedical, Inc. Striate+™ is manufactured by Orthocell Ltd. Striate+™ and SMRT™ are trademarks of Orthocell Ltd. SynMax®, PermaPro® and DEDICAM® are registered trademarks of CAMLOG Biotechnologies GmbH. They may, however, not be registered in all markets.

BioHorizons and MinerOss® A products are approved for sale in the European Union in accordance with pharmaceutical legislation, the Medical Device Directive 93/42/EEC (and where applicable, Regulation 2017/745) and the Human Tissues and Cells Directive 2004/23/EC, respectively. We are registered to ISO 13485:2016, the international quality management system standard for medical devices, which supports and maintains our product licenses with Health Canada and in other markets around the globe. All rights reserved. Not all products shown or described in this literature are available in all countries.

